

pH-Fix 高精度試紙

多年來，未經培訓的用戶和分析專業人員都喜歡使用pH-Fix輕鬆進行pH測試的體驗。
pH-Fix 試紙條中的染料會化學結合到試紙上，
這項專利技術可防止染料滲出，從而防止樣品污染。
試紙採用塑膠片固定色塊可以長時間進行安全的測定。





試紙影片介紹

Test	REF	Measuring range (visual)	Measuring range (instrumental) ¹⁾
0-14	92110	0.1·2·3·4·5·6·7·8·9·10·11·12·13·14	1-13
0.0-6.0	92115	0.0·0.5·1.0·1.5·2.0·2.5·3.0·3.5·4.0·4.5·5.0·5.5·6.0	0.5-6.0
2.0-9.0	92118	2.0·2.5·3.0·3.5·4.0·4.5·5.0·5.5·6.0·6.5·7.0·7.5·8.0·8.5·9.0	2.0-9.0 胃液
4.5-10.0	92120	4.5·5.0·5.5·6.0·6.5·7.0·7.5·8.0·8.5·9.0·9.5·10.0	4.5-10.0 尿液 唾液
6.0-10.0	92122	6.0·6.4·6.7·7.0·7.3·7.6·7.9·8.2·8.4·8.6·8.8·9.1·9.5·10.0	6.0-10.0
7.0-14.0	92125	7.0·7.5·8.0·8.5·9.0·9.5·10.0·10.5·11.0·11.5·12.0·12.5·13.0·13.5·14.0	7.0-13.5
0.3-2.3	92180	0.3·0.7·1.0·1.3·1.6·1.9·2.3	-
1.7-3.8	92190	1.7·2.0·2.3·2.6·2.9·3.2·3.5·3.8	-
3.1-8.3	92135	3.1·3.5·3.9·4.3·4.7·5.1·5.5·5.9·6.3·6.7·7.1·7.5·7.9·8.3	- 唾液
3.6-6.1	92130	3.6·4.1·4.4·4.7·5.0·5.3·5.6·6.1	3.6-6.1 陰道分泌物、透析
4.0-7.0	92137	4.0·4.4·4.7·5.0·5.3·5.5·5.8·6.1·6.5·7.0	- 陰道分泌物
5.1-7.2	92140	5.1·5.4·5.7·6.0·6.3·6.6·6.9·7.2	-
6.0-7.7	92150	6.0·6.4·6.7·7.0·7.3·7.7	6.0-7.7
7.5-9.5	92160	7.5·7.9·8.2·8.4·8.6·8.8·9.1·9.5	-
7.9-9.8	92170	7.9·8.3·8.6·8.9·9.1·9.4·9.8	-
0-14 PT	92111	0.1·2·3·4·5·6·7·8·9·10·11·12·13·14	1-13
3.6-6.1 PT	92131	3.6·4.1·4.4·4.7·5.0·5.3·5.6·6.1	3.6-6.1
4.5-10.0 PT	92121	4.5·5.0·5.5·6.0·6.5·7.0·7.5·8.0·8.5·9.0·9.5·10.0	4.5-10.0

100 test/box

Using pH-Fix test strips

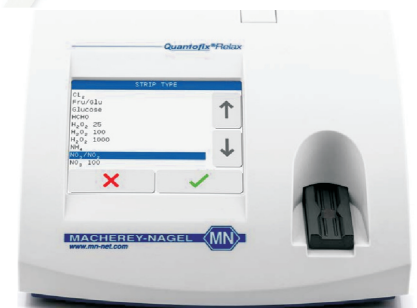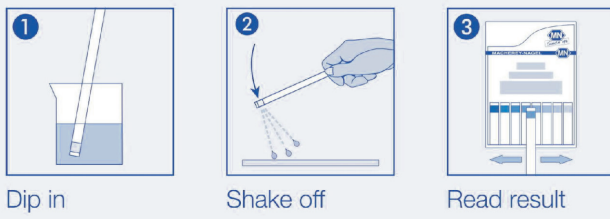

用於評估試紙條的反射儀

Relax 91346

PEHANON直讀式試紙

有色樣品的pH測定。

PEHANON® 試紙在一條試紙上統一了 pH 指示劑和參考色卡。

任何樣品顏色對參考顏色和反應墊都具有相同的效果，

即使在有色溶液中也能獲得純正的 pH 讀數。

頂部色域上方的隱形疏水屏障可防止樣品遷移。

PEHANON® 試紙無需單獨的色卡即可讀取。

生產工人可以使用單條而不是完整的包裝，使產品非常經濟。



試紙影片介紹

Test	REF	Measuring range	Number of Tests
■ pH 1-12	90401	1·2·3·4·5·6·7·8·9·10·11·12	200
■ pH 0-1.8	90411	0·0.3·0.6·0.8·1.0·1.2·1.5·1.8	200
■ pH 1.0-2.8	90412	1.0·1.3·1.6·1.8·2.0·2.2·2.5·2.8	200
■ pH 1.8-3.8	90413	1.8·2.1·2.4·2.7·3.0·3.2·3.5·3.8	200
■ pH 2.8-4.6	90414	2.8·3.1·3.4·3.6·3.8·4.0·4.3·4.6	200
■ pH 3.8-5.5	90415	3.8·4.0·4.2·4.4·4.6·4.9·5.2·5.5	200
■ pH 4.0-9.0	90424	4.0·4.5·5.0·5.5·6.0·6.5·7.0·7.5·8.0·8.5·9.0	200
■ pH 5.2-6.8	90416	5.2·5.5·5.7·5.9·6.1·6.3·6.5·6.8	200
■ pH 6.0-8.1	90417	6.0·6.3·6.6·6.9·7.2·7.5·7.8·8.1	200
■ pH 7.2-8.8	90419	7.2·7.4·7.6·7.8·8.0·8.2·8.5·8.8	200
■ pH 8.0-9.7	90420	8.0·8.2·8.4·8.6·8.8·9.1·9.4·9.7	200
■ pH 9.5-12.0	90421	9.5·10.0·10.5·11.0·11.5·12.0	200
■ pH 10.5-13.0	90422	10.5·11.0·11.5·12.0·12.5·13.0	200
■ pH 12.0-14.0	90423	12.0·12.5·13.0·13.5·14.0	200



通用試紙 Universal and indicator papers

Universal pH試紙 – 許多應用的標準。

對於每個 pH 值，這些紙顯示單一顏色，可以以 0.2~1 pH 單位的間隔與色標匹配。

試紙採用塑料捲筒包裝，可確保長期穩定性並防止許多外部影響。

刻度的顏色經過特殊混合，與試紙的反應顏色完美匹配，

這使得結果的讀取變得容易和準確。



試紙影片介紹



Test	REF	REF refill	Measuring range	Presentation
■ pH 1–11	90201	90202	1·2·3·4·5·6·7·8·9·10·11	Reel of 5 m length and 7 mm width
■ pH 1–11	90203	–	1·2·3·4·5·6·7·8·9·10·11	Booklet with 100 strips 10 x 70 mm
■ pH 1–14	90204	90224	1·2·3·5·6·7·8·9·10·12·14	Reel of 5 m length and 7 mm width
■ pH 0.5–5.5	90205	90225	0.5·1.0·1.5·2.0·2.5·3.0·3.5·4.0·4.5·5.0·5.5	Reel of 5 m length and 7 mm width
■ pH 3.8–5.8	90206	90226	< 3.8·3.8·4.1·4.3·4.5·4.7·4.9·5.2·5.5·5.8·> 5.8	Reel of 5 m length and 7 mm width
■ pH 4.0–7.0	90207	90227	4.0·4.3·4.6·4.9·5.2·5.5·5.8·6.1·6.4·6.7·7.0	Reel of 5 m length and 7 mm width
■ pH 5.4–7.0	90208	90228	< 5.4·5.4·5.7·6.0·6.2·6.4·6.7·7.0·> 7.0	Reel of 5 m length and 7 mm width
■ pH 5.5–9.0	90209	90229	5.5·6.0·6.5·7.0·7.5·8.0·8.5·9.0	Reel of 5 m length and 7 mm width
■ pH 6.4–8.0	90210	90230	< 6.4·6.4·6.6·6.8·7.0·7.2·7.4·7.6·7.8·8.0·> 8.0	Reel of 5 m length and 7 mm width
■ pH 7.2–9.7	90211	90231	< 7.2·7.2·7.5·7.8·8.1·8.4·8.7·9.0·9.3·9.7·> 9.7	Reel of 5 m length and 7 mm width
■ pH 8.0–10.0	90212	90232	8.0·8.2·8.4·8.7·9.0·9.2·9.6·10.0	Reel of 5 m length and 7 mm width
■ pH 9.0–13.0	90213	90233	9.0·9.5·10.0·10.5·11.0·11.5·12.0·12.5·13.0	Reel of 5 m length and 7 mm width
■ pH 12.0–14.0	90214	90234	12.0·12.5·13.0·13.5·14.0	Reel of 5 m length and 7 mm width

新品上市測油脂酸價



遊離脂肪酸試紙(油脂老化)

產品編號：REF 91365

量測範圍：0.25~3% FFA

範圍區間：0.25-0.5-1-2-3%FFA

鋁製容器內含50片試片